
Aukey International Limited
Web: www.naipocare.com
E-mail: wholesale@naipocare.com
Add.: Room 1202 Hong Commercial Building
145 Hennesy Road Wanchai Hong Kong
HQ: No. 102, Building P09, China South City, No. 1 Huanan Avenue,
Pinghu Street, Longgang District, Shenzhen, China

IP44
Made in China

User Manual

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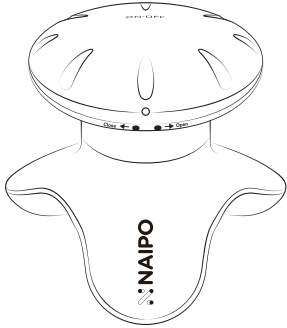Model: MGPC-101MM

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Introduction

Thank you for purchasing the Naipo Mini Massager. Please read this user manual carefully and keep it in a safe place for future reference.

Included in the Package

1 x Vibration Massager
3 x AAA Batteries
1 x User Manual

Product Specifications


Model Name	Naipo Mini Massager
Weight	≈0.26 kg
Size	11.2 x 10.3 x 9.3 cm
Material	ABS, Resin gum
Power Input	DC 4.5 V 3 x AAA Batteries
Power Consumption	1.3 VA
Massage Intensity	Low
Water resistant	Protection rating IP44

- EN -

Product Operation Guide

Batteries Installation

1. Open the cover of the massager and install 3 AAA batteries according to the indicated direction in this diagram:

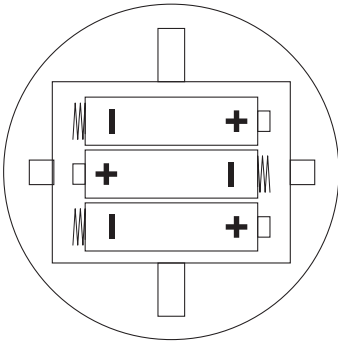2. Confirm that the anode and cathode orientation of the batteries is correct and make sure that the battery cover is properly closed.

- EN -

Operation

- Place the 3 massage nodes on the part of your body that you want to massage and press downwards the top of the massager.
- The massager will vibrate as long as the top is pressed with enough pressure.

Troubleshooting

Problem	Possible cause	Solution
The device doesn't work.	Low batteries or Wrong direction of the anode and cathode of the batteries	Change the batteries or check if the direction of the anode and cathode is correct
The vibration is abnormally weak	Low batteries	Change the batteries

- EN -

Maintenance

Moisten a soft cloth with alcohol and wipe the surface of the massager. Please keep your massager on a safe place after use. Do not place any objects on top of your massager. Do not disassemble.

Transportation and Storage

Ambient temperature: 15°C - 35°C
Maximum relative humidity: 75%
Ventilated room and no corrosive gas
Atmosphere pressure: 860 mbar - 1060 mbar
The massager is lightly shockproof., suitable for airplane, railway, road, and ship transportation.
Do not expose to immersion or collision.

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