
### User Manual
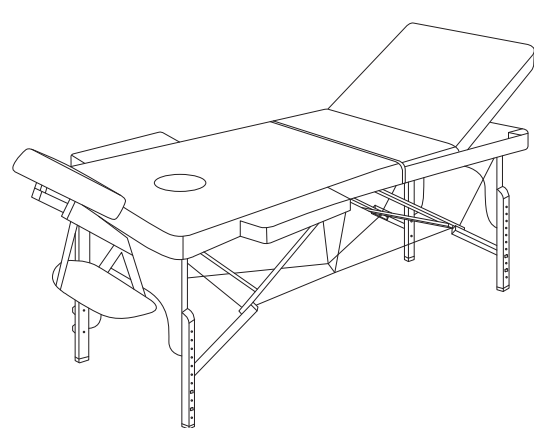
### Three-fold Portable Massage Table with Wooden Feet

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#### Introduction

Thank you for purchasing the Naipo Three-fold Portable Massage Table with Wooden Feet. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.


#### Included

- 1 x Massage table
- 1 x Face hole plug
- 1 x Face cradle cushion
- 1 x Face cradle
- 1 x Armrest
- 2 x Side armrest
- 1 x Carrying bag
- 1 x User manual

#### Product Specifications

Model	MGBC-302H
Name	Naipo Three-fold Portable Massage Table with Wooden Feet
Dimensions	185*70*(63-88)cm (adjustable height)
Maximum working weight capacity	270kg

#### Using Guide

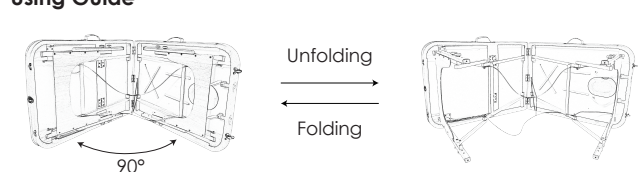
- 01 -

#### User tips:

- When unfolding, lay the folded product on its side, unfold it until the angle is 90°, and then unfold the wooden feet.
- When folding, lay the unfolded product on its side, fold it until the angle is 90°, and then fold the wooden feet.
- Please make sure that the wire ropes do not wind around the wooden feet when unfolding or folding the product.


#### Installation of Face Cradle Cushion

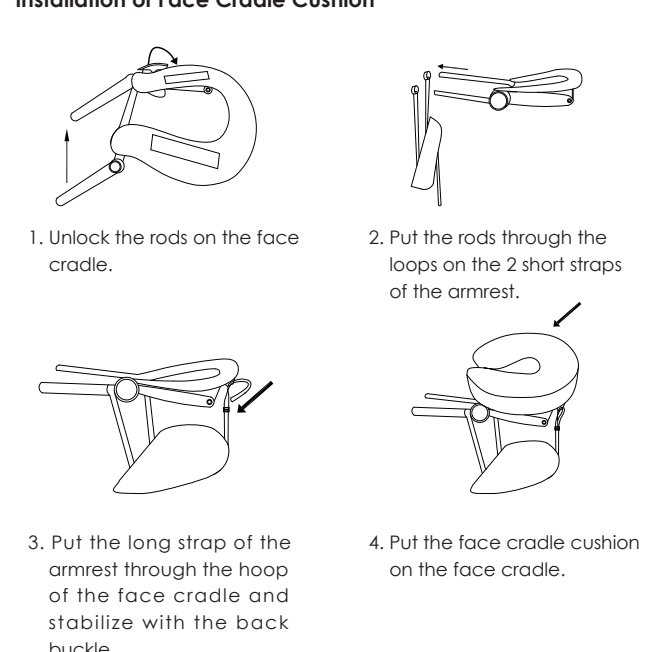- 02 -

#### Height Adjustment


**Instructions:** Rotate the 2 knobs counterclockwise on the wooden foot until they fall down. Take down the wooden foot and put it on the desired height with the screws through the corresponding holes. Then put the knobs back and stabilize by rotating them clockwise.
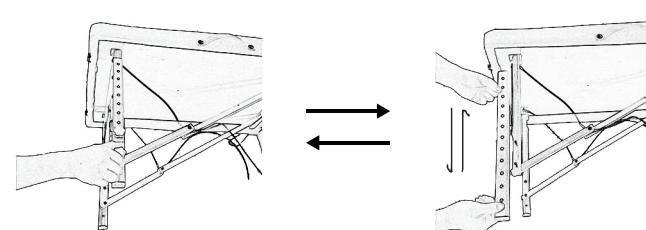
#### Adjustable Backrest

**Lifting:** Gently lift the backrest with 2 hands, and the lifter would respond accordingly, which allows the selection of desired height.

**Lowering:** Lift the backrest to the maximum height and then put it down.

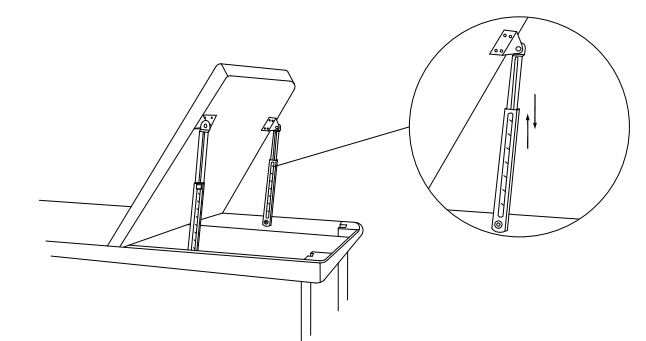- 03 -

#### Product Maintenance

- The product is designed for massage, not for intensive activities.
- Please wipe the product surface with neat sponge or mildly damp cloth for cleaning.
- Please do not clean the product surface with any corrosive detergent.
- Please fold and store the product if it's not to be used for long.

#### Warranty & Customer Support

If you need help regarding product issues, after-sales service, or other relevant problems, please contact us via the email on the cover.

\* Note: Products purchased directly from Naipo will enjoy our after-sales service. If the purchases are made from other selling teams, we recommend you contact them for warranty issues.

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