Foot and Calf 2 in 1 Massage with Air Pressure & Rolling Massage

Model: MGF-3600

Shenzhen Pinlong E-Commerce Co., Ltd.
Web: www.naipocare.com
E-mail: support@naipocare.com
Add: Huanan City Electronic Trading Building P07 Room 120, Longgang District, Shenzhen, China

Introduction

Thank you for purchasing the Naipo Foot and Calf 2 in 1 Massage with Air Pressure & Rolling Massage. Please read this manual carefully to ensure that it is used safely and effectively. If you need any assistance, please contact our support team with your product model number and Amazon order number.

- Included**
 - 1 x Foot & Calf Massager
 - 1 x Adapter (Please use only with the designated adapter)
 - 1 x User Manual

Product Specifications	
Model	MGF-3600
Name	Naipo Foot and Calf 2 in 1 Massage with Air Pressure & Rolling Massage
Input Voltage	DC 24V
Power consumption	55W
Auto-off time	15 minutes

Important Safeguards

- Close supervision is necessary if the product is to be used by those who are disabled, who suffer from sensory or mental incapacity, and who lack the relevant experience and knowledge of the product (children included).
- Keep the product away from children.
- The surface of the product might be heated due to the heating device inside. Therefore, people who are insensitive to heat should be aware of this.
- To ensure safety, if the power cable is damaged, consult with the manufacturer, designated maintenance institutions, or professional maintenance teams to replace the cable.
- Keep the product away from sharp objects. Do not use the product if fabric covers are damaged so that foot holes do not occur. Consult with the manufacturer, designated maintenance institutions or professional maintenance teams to repair the product before using.
- Please disconnect from the power supply when the product is not in use or before it is to be cleaned.
- Do not cover the product while it's in use, otherwise the motor will get overheated, which might cause malfunction or damage to the product.
- Stop using the product immediately if uncomfortable feelings occur during use.
- If the product malfunctions, please consult with the manufacturer, designated maintenance institutions or professional maintenance teams to repair the product before using.
- Use only spare parts recommended by the manufacturer.
- Do not drop or insert objects into any opening of the product.
- Do not operate the product with wet hands.
- Do not contact the control panels with water or other liquids.
- Turn off the product before disconnecting from the power supply.
- The product is not intended for wet diagnosis or self-treatment for conditions that should be managed by a qualified health care provider.
- It's recommended to use the product in an environment with temperature between 5°C (41°F) - 40°C (104°F), and please do not use it in a humid environment.
- Do not use the product outdoors or expose it to direct sunlight.

Warnings

Do not move the product from low-temperature to high-temperature, for the water vaporization might damage the product. After using the product, it's recommended that it be put in a place at room temperature for 1 hour.


- People listed below should consult a doctor before using the product:**
- People who are using a medical treatment or people who are feeling uncomfortable.
 - People who are suffering from hypertension, heart disease, central arteriosclerosis, or cerebral thrombosis.
 - People who have a malignant tumor.
 - Women who are pregnant or women who are in menopause.
 - People who are suffering from osteoporosis.
 - People who are misdiagnosed with gonorrhea or other medical apparatus.
 - People who are suffering from an acute disease.

Q&A

If the product fails to work after the measures below have been taken, please send the product to designated service center or designated technical staff for maintenance or repair. Note: It's recommended to use the massager for 1-2 times per day, with each time less than 15 minutes.


Problem	Cause	Remedy
The product doesn't work.	The product is not connected to the power supply.	Please connect the product to the power supply or make sure that the connection is good.
	The power button has not been pressed.	Please press the power button.
The product stops working while in use.	The product is not connected to the power supply tightly.	Please connect the product to the power supply tightly.
	The product has been in consistent use for a long period of time, which might activate overheat protection.	Please press the power button again.
The product has been used for long periods of time, which might activate overheat protection.	The pressure is too high for the massage nodes.	Adjust the load by relaxing feet.

Product Diagram

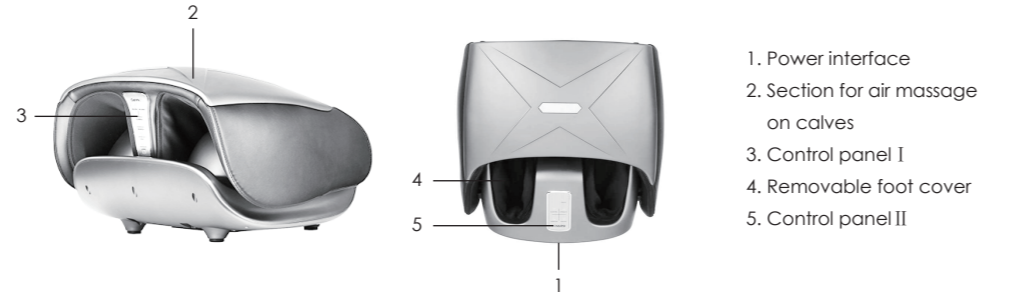
Control Panel I

Press **U** to turn on or turn off the massager. When it's turned on, rolling massage on feet will be activated by default.
Press **1** to control air compression among 3 intensity levels.
Press **2** to control vibration function on foot among 3 intensity levels or turn off the function.
Press **3** to activate or turn off the heating function.
Press **4** to activate or turn off leg mode 1 / leg mode 2.
Press **5** to activate or turn off foot mode 1 / foot mode 2.
Under leg mode, long press **4** to stop rolling massage on feet, while the air compression on sides of calves, calves, and feet, while the air compression on sides of calves, calves, and feet remains. The intensity could still be controlled. Press **3** again to start to normal working mode.
Under foot mode, long press **2** to stop rolling massage on feet, while the air compression on sides of calves, calves, and feet remains. The intensity could still be controlled. Press **3** again to start to normal working mode.
Under leg mode, long press **1** to stop air compression while the rolling massage on feet remains.
Under foot mode, long press **1** to stop air compression while the rolling massage on feet remains.

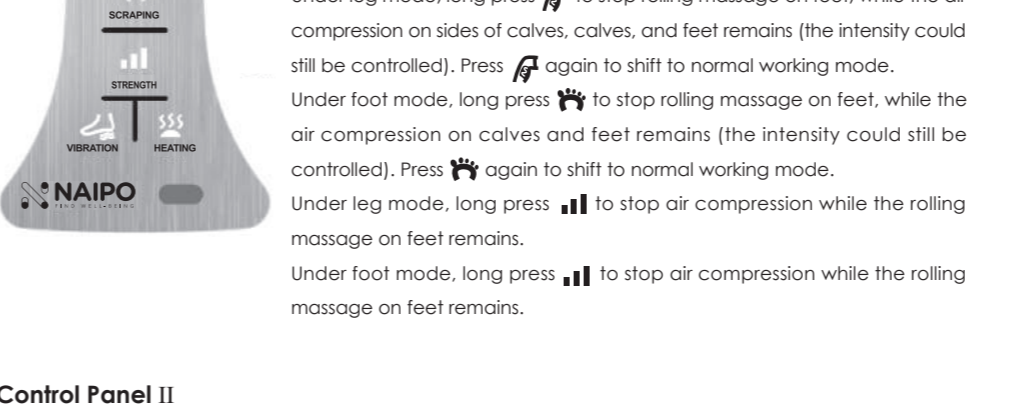
Control Panel II

Press **U** to turn on or turn off the section for foot massage. When it's turned on, rolling massage on feet, air compression on sides of calves, calves, and feet, and vibration on feet will be activated by default.
Press **1** to control air compression among 3 intensity levels.
Press **2** to control vibration function on foot among 3 intensity levels or turn off the function.
Press **3** to activate or turn off the heating function.
Press **4** to activate or turn off leg mode 1 / foot mode 2.
Long press **2** to stop rolling massage on feet while the air compression on feet remains. The intensity could still be controlled. Press **3** again to start to normal working mode.
Long press **1** to stop air compression on feet while the rolling massage on feet remains. Press **3** again to activate air compression on feet again.

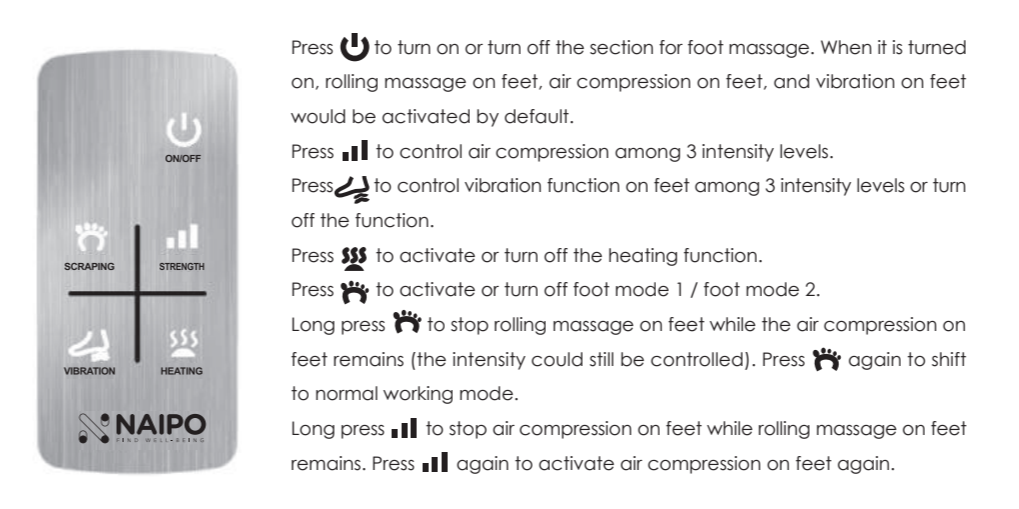Product Operation Guide

- Place the product flat on a clean floor. Make sure that the foot covers are not placed without any damage.
- Connect the product to the power supply with the adapter and the product will start to standby mode with the power indicator flashing. Sit on a chair and put legs inside the foot covers (do not stand while using the product).
- When faded, the product can massage on feet. When unfolded, the product can massage on feet and calves. Choose the power you need with different massage modes.
- Control panel I only works when the product is unfolded, while Control panel II only works when the product is folded.
- Press turn off the product if it is to be unfolded or folded.
Note: When folding, please turn off the product and leave it alone for 1 minute until the air inside the air bag of the calf action has been completely released.
- While in use, the product will shift into standby mode if the power button is pressed or the auto-off function has been activated.

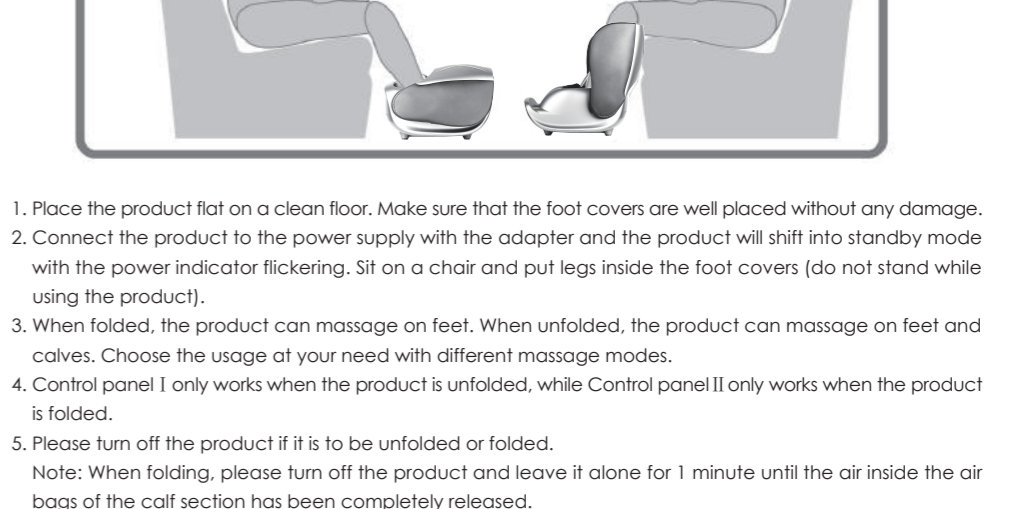Massage Areas

This product is designed for massage on feet and calves.

Warranty and Customer Support

For warranty, assistance, support, or warranty claims, please contact us at the address below that corresponds with your region. Please include your Amazon order number and product model number.
Address: Huanan City Electronic Trading Building P07 Room 120, Longgang District, Shenzhen, China.
Website: www.naipocare.com
Email: support@naipocare.com

*Note: Products purchased directly from Naipo cannot be returned for warranty issues. If the purchase is made from other selling teams, we recommend you contact them for after-sales service.

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